2 The synergy of PDRN's regenerative and repairing action with microneedling technology stimulates the skin cells' natural renewal mechanisms and collagen production, while also enabling immediate and faster penetration of PDRN and other serum active ingredients.


3 The synergy of PDRN with the following targeted active ingredients to enhance and improve key vital functions of the skin and its physiology:

- Copper Tripeptide enhances collagen production and, in combination with the botox-like hexapeptide Argireline, improves skin firmness and elasticity.
- B-complex Vitamins (B3, B5, B12) and stabilized Vitamin C (AA2G) provide strong antioxidant and protective action, shielding the skin from free radicals and environmental pollutants, improving discolorations, and offering radiance.
- Hyaluronic acid (MMW), along with a specialized complex of moisturizing ingredients, maintains and enhances skin moisture levels.
- Centella Asiatica and Vitamin B12 act soothingly and regeneratively, preventing and alleviating potential irritations.

Application Instructions

1 Cleansing: Clean your skin thoroughly with an appropriate cleanser according to your skin type.


Press the head down once to release the appropriate amount of serum.

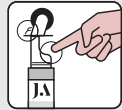
2 Treatment Areas: Divide the face into 6 areas, applying the product to each area, focusing more on areas with greater skin damage (e.g. wrinkles, sagging skin, discolorations).

Optionally, the product can also be applied on the neck, dividing it into two areas (left and right), avoiding the thyroid area.

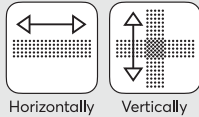
Horizontally

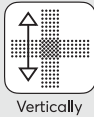Vertically

3 Application: After disinfecting the needles, apply the product on each area as follows:

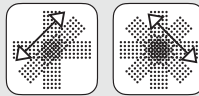Diagonally

4 Calming Bliss: Apply the Calming Bliss for immediate skin soothing. In case of intense irritation, apply frequently for the first 24 hours.

Place the roller with the needles on each area, immobilize the skin with one hand, and move the roller with the other hand, without applying much pressure, along the skin: horizontally, vertically, and diagonally (4 times per direction), lifting the device each time.

5 Sun Protection: Throughout the application of the product, the use of daily high sun protection is recommended.


The product is recommended for 10 applications on the face, one application at least every 7 days. Alternatively, it can be used both on the face and neck for 8 applications.


Cleaning & Disinfecting the device

1 After use, clean the head under running water with the special brush.

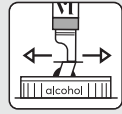2 Soak the sponge of the auxiliary disinfecting jar in alcohol and disinfect the needles (before and after use).

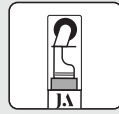3 Allow the head to dry well and put the cap on.

Contraindications

Do not use it in the following cases: irritated skin, skin infections and malignancies, acne, moles, wounds, pregnancy, breastfeeding.

Caution

For external use only, exclusively by one individual.

After microneedling, it is recommended to avoid exposure to ultraviolet (UV) radiation, and the use of saunas and steam baths (hammam).

Side Effects

Side effects have not been documented with the proper use of the PDRN Rejuvenation Microneedling.

Do not hesitate to contact a professional or the company regarding the instruction leaflet or the product.

Storage

Store the product in a safe place, out of reach of children. Store in a dry environment, away from sunlight.

Home Daily Regimen

For your daily skincare regimen (morning and night), use the PDRN Rejuvenation Serum, and the appropriate cream based on according to your skin's needs, DNA-RNA Rejuvenation Cream or PDRN Rejuvenation Cream.

PDRN Rejuvenation Serum
Skin rejuvenation face serum with PDRN polynucleotide, enriched with peptides, vitamins and hyaluronic acid. For all skin types.

Indications: Wrinkles, expression lines, sagging skin, scars, discolorations, dull skin

DNA-RNA Rejuvenation Cream
Skin rejuvenation face cream with a complex of hydrolyzed DNA-RNA nucleotides, enriched with hyaluronic acid, peptides and vitamins. For all skin types.

Indications: Fine lines, expression lines, scars, dull skin

PDRN Rejuvenation Cream
Skin rejuvenation face cream with PDRN polynucleotide, enriched with peptides, ceramides, vitamins and hyaluronic acid. For normal, dry and mature skin.

Indications: Wrinkles, expression lines, sagging skin, scars, discolorations, dull skin.

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
EN

PDRN Rejuvenation Microneedling

Advanced skin rejuvenation system for the face combining microneedling technology with a PDRN (PolyDeoxyRiboNucleotide) polynucleotide serum enriched with peptides, vitamins, and hyaluronic acid. For all skin types.

The controlled microtrauma caused during the application stimulates skin's natural renewal mechanisms with immediate and faster penetration of the active ingredients. The combined action of microneedling technology in synergy with the ingredients of the serum regenerates the skin, enhancing its elasticity and improving texture and skin tone.


The box contains:

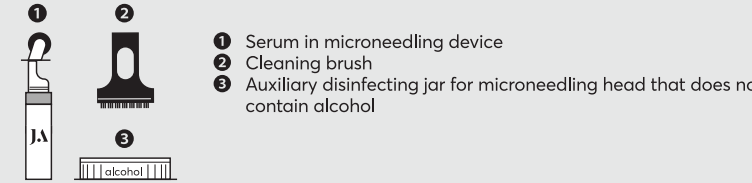
Indications

Wrinkles, expression lines, sagging skin, scars, discolorations, dull skin

Characteristics

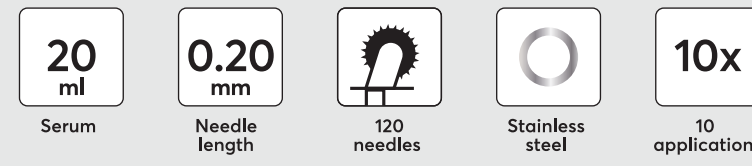Mechanism of Action

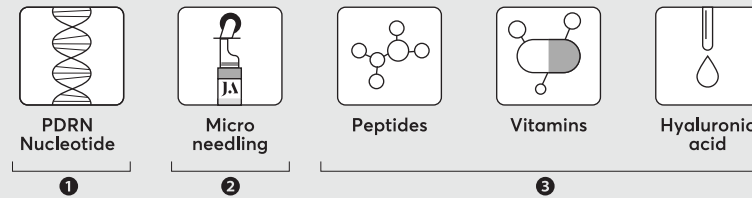The action of PDRN Rejuvenation Microneedling is based on:

- 1 The multifaceted mechanism of PDRN polynucleotide in a biological and mechanical way:
 - It binds to receptors on epidermal cells and fibroblasts, promoting cell proliferation and collagen and elastin production, reducing inflammation and oxidative stress, ensuring regular and complete cellular regeneration.
 - It significantly contributes to the skin's cells DNA "self-repair" mechanism as a source of free nucleotides.
 - Its polymeric organic structure forms a film on the skin's surface, offering an immediate tightening and "lifting" effect.